


<b>Product:</b>	<b>MAMMOMAT Novation DR (DROC)</b>
<b>Title:</b>	<b>Specific Hardcopy Camera Information for Kodak DryView 8610 / MIM 200</b>
Affected systems: MAMMOMAT Novation DR with Image Processing, V1.0 and Hardcopy Cameras	
<b>Replaces:</b>	n.a.

## 1 General Information

### 1.1 Block Diagram

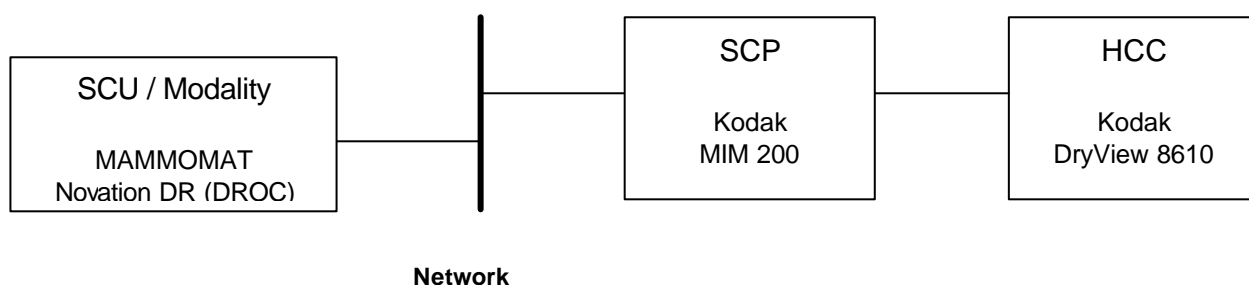

### 1.2 Available Film Size

Application			Film Type	Medium Type
14"x17"	Portrait	n.a.	n.a.	n.a.
	Landscape	n.a.		
14"x14"	Portrait	n.a.	n.a.	n.a.
	Landscape	n.a.		
11"x14"	Portrait	n.a.	n.a.	n.a.
	Landscape	n.a.		
10"x12"	Portrait	n.a.	n.a.	n.a.
	Landscape	n.a.		
8"x10"	Portrait	yes	DVM	Blue-based
	Landscape	yes		
others		n.a.	n.a.	n.a.

yes - Film size is available and released  
no - Film size is available but not released  
n.a. - Film size is not available

#### 1.2.1 Supported Format

Supported page format is always 1x1

# MAMMOMAT Novation DR (DROC) Camera Installation

## Specific Kodak DryView 8610 / MIM 200 Information

### 1.3 Restrictions

For HCC-independent restrictions, see also the general restrictions in the MAMMOMAT Novation DR document.

The linear LUT of max. density 3.5 is out of tolerance (in the basic qualification), but for mammography images the expected reference LUT is reached.

## 2 Modality SW Configurations / Parameters

Refer to: MAMMOMAT Novation DR/CB-DOC/ Installation/Installations instructions.

file name	content	remarks
dryview8610.cfg	<p><b>Configuration DROC</b> Software Configuration Edit ./installed/raw/vdevdryview8610_001 Edit ./installed/raw/vdevdryview8610_001's configuration file: ./linx_mp/resources/installed/raw/vdevdryview8610_001/dryview8610.cfg</p> <p>-----</p> <p><b>Group: Server</b> factory object name = RMS://siemens07/installed/services/jss_server:PJF_OREF spooler name = DicomPrintSpooler priority = normal max retries = 20</p> <p>-----</p> <p><b>Group: Output Device</b> device type = DRYVIEW</p> <p>-----</p> <p><b>Group: Print Job Description</b> device name = DV8610 device address = 146.254.124.43:5040 request annotation SOP class = false request print job SOP class = false print job poll period = 15 print job timeout = 300 printer poll period = 15 association timeout = 0 collate = false custom lut file = photometric interpretation = automatic</p> <p>-----</p> <p><b>Group: Film Session</b> copies = 1 medium type = BLUE FILM film destination = PROCESSOR film session label = DROC film memory allocation = 0 print priority = MED</p> <p>-----</p> <p><b>Group: Film Box</b></p>	<p>This parameter can be changed within Netscape Service Login. Only imported parts of the file are shown.</p>

# MAMMOMAT Novation DR (DROC) Camera Installation

## Specific Kodak DryView 8610 / MIM 200 Information

	<p>             available image display formats = PORTRAIT\1,1              image display format = PORTRAIT\1,1              annotation display format id = 6              number of annotation boxes = 6              annotation boxes per text line = 3              max annotation box text string size = 64              annotation string =              max text line size = 189              number of image text lines = 5              image text position = left              image text overlap = overwrite              image text reorient = reorient_0  <b>film size id = 8INx10IN</b>  <b>film magnification type = REPLICATE</b>              film smoothing type = NONE              border density = BLACK              empty image density = BLACK              full pixel range = true              min density = 20  <b>max density = 350</b>  <b>configuration info = LUT=0,5</b>              trim = NO           </p> <p>-----</p> <p> <b>Group: Image Box</b>              polarity = NORMAL              available image magnification types = REPLICATE BILINEAR              CUBIC  <b>image magnification type = REPLICATE</b>              image smoothing type = NONE  <b>true size capable = true</b>  <b>enable selective true size = true</b>              pixel size = 70              pad = false              image text line 1 = &lt;Off&gt;              image text line 2 = &lt;Off&gt;              image text line 3 = &lt;Off&gt;              image text line 4 = &lt;Off&gt;              image text line 5 = &lt;Off&gt;           </p> <p>-----</p> <p> <b>Group: PORTRAIT\1,1</b>              width = 2781              height = 3430           </p>	
--	---	--

Note:	After changing Film Size ID in Netscape, you have to press the Update button and restart the Sun!!!
-------	---

# MAMMOMAT Novation DR (DROC) Camera Installation

## Specific Kodak DryView 8610 / MIM 200 Information

### 3 Hardcopy Camera, Configuration

#### 3.1 Identification

<b>Device Type:</b>	<b>DICOM Print Server(SCP)</b>	<b>Printer</b>
<b>Manufacturer / Device Name:</b>	Kodak / PACS Link MIM 200	Kodak / DryView 8610
<b>HW Revision:</b>	IBM – PC, xSeries 200	n.a.
<b>SW Revision:</b>	V 6.1	V 2.4
<b>Image Memory:</b>	256 MByte	256 MByte
<b>Hard disk:</b>	4 GByte	n.a.
<b>Number of Operating Panels:</b>	1	1 Built In
<b>No./Type of Input:</b>	1 Ethernet / RJ45	1 Control / 1 Data
<b>No./Type of Output:</b>	1 Copcil Board	n.a.
<b>Physical Network Interface:</b>	10/100 Base T/Tx	n.a.
<b>Network Protocol:</b>	TCP/IP	n.a.
<b>DICOM Conformance Statement:</b>	Kodak PACS Link MIM 200 / Software Version 6.1 / October 17, 2003 / Rev. G	n.a.

#### 3.2 SCP / Printer Settings

File name	Content	Remarks
Not available	<b>MIM200 Setup:</b> <ul style="list-style-type: none"> <li>■ on Configure, Network:</li> <li>■ Node Name=MIM, IPAddr=accordingly</li> <li>■ Subnet Mask=accordingly, Gateway Addr=accordingly</li> <li>■ on Configure, Input, Network, Add:</li> <li>■ Protocol=DICOM 3.0, Log.Name=DICOM, Type/Com=Keypad/COM9</li> <li>■ on Configure, Destination, Add, Print, Direct, Qualified, 8610, Set:</li> <li>■ General: Log.Name=8610</li> <li>■ General: Database Name=8610, Model=8610, Memory=32</li> <li>■ Media Info: Media=8x10in, 5025x6200, 6020, Blue</li> <li>■ Output Params: Cmd Set=952, Baud=1200, Data=8, Stop=1, Parity=Even</li> <li>■ Output Params: Pixel=12, Header/Line, Transfer=2000 KHz, Acq=60</li> <li>■ Supported Item: Collation=On, Text=On, Tonescalning=On, Rotation=On</li> <li>■ Interpolation: NoMag=On, Bilin=On, Repl=On, Cubic Spline=On</li> <li>■ Format Types: Custom=On, Slide=On, Row Symmetric=On</li> <li>■ Format Types: 1, 2, 4, 6, 8, 9, 12, 15, 16, 20, 24, 30, 35=On</li> <li>■ Formatting: Aspect Ratio Mix on Page=On, Image Size Mix in Std Fmt=On</li> <li>■ Formatting: DICOM Annotation Box=On, Req.Size Scale=On</li> <li>■ Formatting: Req.Size Mix in Row=On, Trim=On</li> <li>■ Formatting: Aspect Ratio Mix in Row=On, Image Size Mix in Row=On</li> <li>■ Formatting: Req.Image Size=On, Req.Size Mix on Page=On</li> <li>■ Formatting: True Landscape=On</li> </ul>	Provided by Kodak

# MAMMOMAT Novation DR (DROC) Camera Installation

## Specific Kodak DryView 8610 / MIM 200 Information

	<ul style="list-style-type: none"> <li>■ Formatting: Max Img per Band=20, Max Band per Page=20, Max Img Col=0</li> <li>■ Formatting: Max Img Row=0, Min Horiz Sep=2, Min Vert Sep=2</li> <li>■ Formatting: Trim Width Enabled=2, Trim Width Disabled=0</li> <li>■ Formatting: Pixel Pitch=25.8065</li> <li>■ Image Processing: Grayscale=On, Minify=On, Max Image Dmax=3.5</li> <li>■ M952 Usage: Min H/V Mag=0.26, Max H/V Mag=16, Min Dmax 952=2</li> <li>■ on Configure, SCP Services, select the specific modality, Modify:</li> <li>■ PLUT=N, ToneScaling=PerceptionLUT (Always=Y)</li> <li>■ <b>Contrast (dry)=5 (Always=N)</b></li> <li>■ <b>Modality Group=1 (which corresponds to TFT/ULUT=upmsv3kp.w8u)</b></li> <li>■ Conv.Group=1, Gamma Table=Linear in Density</li> <li>■ Image Dmin=0.00 (Always=N), Image Dmax=3.50 (Always=N)</li> <li>■ Smoothing Type=1 (Always=N), Image Polarity=Normal (Always=N)</li> <li>■ NeventReports=Off (Always=Y), Warning Status=On</li> <li>■ Print Date=Off, Print Time=Off, Page Number=Off, Hosp.Name=Off</li> <li>■ Modality ID=Off, Annotation Override=Off (Always=N)</li> <li>■ Minify/Scale=On, Crop=Off, Source ID=droc</li> <li>■ Print Priority=Normal (Always=N), Border=3.99 (Always=N)</li> <li>■ Film Base=Blue Film (Always=N), Film Size=8INx10IN (Always=N)</li> <li>■ Negotiate Color Print=Off</li> </ul> <p><b>8610 Setup:</b> on 8610 Operator Menu into Test Menu, Software Version, verify:</p> <ul style="list-style-type: none"> <li>■ Imager=2.4, Contrast=6.P.2, Conv.Kernel=M.1.0</li> <li>■ (if 8600: Imager=2.5, Contrast=6.G.7, Conv.Kernel=1.0.5)</li> </ul> <ul style="list-style-type: none"> <li>■ on 8610 Operator Menu into Configure User set: <ul style="list-style-type: none"> <li>■ Density=16, Contrast=6, Sharp/Smooth=1, Film Layout: 1:1P</li> <li>■ Gamma=Xmittance, Image Size=Scale To Fit, Detector Size=don't care</li> </ul> </li> <li>■ on 8610 Service Menu into Configure User set: <ul style="list-style-type: none"> <li>■ Ratio=1.00, Serial Control</li> <li>Intf.=1200,8,1,Even,Comma,3M, Memory=127</li> <li>■ Digital Interface: Parity=Even, Modality Group=1, Smoothing Group=1</li> </ul> </li> </ul>	
--	---	--

### 3.3 Error and Warning Messages

n.a.

## 4 Image Quality and Optical Density

# MAMMOMAT Novation DR (DROC) Camera Installation

## Specific Kodak DryView 8610 / MIM 200 Information

---

### 4.1 Reference Film and LUT

Released Film Type	DVM	BLUE FILM
--------------------	-----	-----------

Refer to: MAMMOMAT Novation DR/CB-DOC/ Installation/Installations instructions.

### 4.2 Real Size Scaling

n.a.

## 5 Abbreviations

AE Title	Application Entity Title
DCS	DICOM Conformance Statement
DICOM	Digital Imaging and Communication in Medicine
HW / SW	Hardware / Software
HCC	Hardcopy Camera
LUT	Look Up Table
OEM	Original Equipment Manufacturer
SCP	Service Class Provider
SCU	Service Class User
TCP/IP	Transfer Control Protocol / Internet Protocol
PLA	Pacsnet Logical Address

## 6 Changes to Previous Version

n.a.

© Siemens AG 2003

The reproduction, transmission or use of this document or its contents is not permitted without express written authority. Offenders will be liable for damages. All rights, including rights created by patent grant or registration of a utility model or design, are reserved.

### Document revision level

The document corresponds to the version/revision level effective at the time of system delivery. Revisions to hardcopy documentation are not automatically distributed.

Please contact your local Siemens office to order current revision levels.

### Disclaimer

The installation and service of equipment described herein is to be performed by qualified personnel who are employed by Siemens or one of its affiliates or who are otherwise authorized by Siemens or one of its affiliates to provide such services.

Assemblers and other persons who are not employed by or otherwise directly affiliated with or authorized by Siemens or one of its affiliates are directed to contact one of the local offices of Siemens or one of its affiliates before attempting installation or service procedures.